### **Original Article**

# Detection of Herpes Simplex Virus Antibodies Using the Whole Virus and Recombinant gD

Meshkat Z<sup>1†</sup>, Soleimanjahi H<sup>1\*</sup>, Roostaee MH<sup>1</sup>, Bamdad T<sup>1</sup>, Meshkat M<sup>2</sup>

- 1. Virology Department, Faculty of Medical Sciences, Tarbiat Modares University. Tehran, Iran.
- 2. Islamic Azad University, Mashhad Branch, Mashhad, Iran.

#### **Abstract**

**Background and Aims:** Herpes simplex virus type 1 (HSV1) remains a potentially serious health problem world wide. All infected people, including asymptomatic ones, are potential sources for virus transmission. Virus envelope contains at least 13 glycoproteins, which glycoprotein D is the major target of immune responses. The aim of this study was development of a specific method that is a more rapid, sensitive and specific test compared to the virus neutralization test which is applied as gold standard test.

**Methods:** In this study, the Western blot technique using crude HSV1 whole particle and baculovirus derived glycoprotein D of HSV1 as antigens was set. Human sera were analyzed by virus neutralization test and then serum samples with reciprocal virus neutralization antibody titers of 32, 64, and 128 were taken to be analyzed by Western blotting.

**Results:** It was shown that there was a very good correlation between results obtained from virus neutralization antibody titers and those of Western blotting. Western blotting using recombinant glycoprotein D of HSV1 as an antigen showed positive results similar to the whole HSV1 antigen.

**Conclusion:** Our study showed the Western blotting using recombinant glycoprotein D can replace virus neutralization test in diagnosis of HSV1 infections.

**Keywords:** Glycoprotein D; Herpes simplex virus type1; Recombinant gD; Virus neutralization test

### Introduction

erpes simplex virus type 1 (HSV1) is an enveloped, double-stranded DNA virus and is one of the most important pathogen of human. It causes a broad range of infections from asymptomatic to life threatening encephalitis. After primary infection, the virus becomes latent in the

\*Corresponding author: Hoorieh Soleimanjahi, Virology Dept., Faculty of Medical Sciences, Tarbiat Modares University, Tehran, Iran. Tel: (+98)21-82 88 35 61; Fax: (+98)21-88 01 30 30; Email: soleim h@modares.ac.ir.

† Microbiology and Virology Research Center and Women Health Research Center, Mashhad University of Medical Sciences, Mashhad, Iran. trigeminal ganglia and persists for the life time (1-4).

The virus genome encodes more than 30 proteins, of which 13 known antigenically distinct glycoproteins are located in the envelope of the virions (5-7).These glycoproteins are the primary inducers of humoral immune response in HSV infection and are also important targets of cell-mediated immune responses (8, 9). Glycoprotein D (gD) is involved in virus infection and pathogenesis (10-14). It is also a suitable candidate for investigation because gD is a confirmed target of neutralizing antibodies, antibody dependent

cell mediated cytotoxicity (ADCC), and virusspecific cellular responses mediated by both CD4+ and CD8+ T lymphocytes (15).

In this study, we tried to run Western blotting using HSV1 whole particles and baculo virus derived HSV-gD1 as antigens and compared the results with virus neutralization test (VNT).

### **Methods**

#### Cell culture

HeLa cells were cultured in Dulbecco's modified Eagle medium (DMEM) supplemented with 5% heat inactivated FCS, 100~IU/ml penicillin and  $100~\mu\text{g/ml}$  streptomycin at  $37^{\circ}\text{C}$  with an atmosphere of 5% CO2.

### Virus propagation

Herpes simplex virus (HSV) was isolated from labial vesicles of a patient and identified as HSV1 using monoclonal antibody (15). HeLa cells were infected at a multiplicity of infection (moi) of 1 pfu/cell, after one hour adsorption, unabsorbed virus suspension was discarded and the cells were washed two times by sterile phosphate buffer saline (pH=7). Serum free DMEM containing 100 IU/ml penicillin and 100 μg/ml streptomycin was added to the cells. After appearance of cytopathic effects (CPE) in almost 70% of the cells, virus was harvested by freezing and thawing. Cell debris was removed by centrifugation at 8000g for 30 min at 4°C and the supernatant was used as whole antigen in Western blotting.

### Virus neutralization test (VNT)

One hundred serum samples were collected from volunteers with or without history of herpetic disease and used for detection of neutralizing antibodies. Briefly, two fold dilutions of heat inactivated of each serum ranging from 1/2 to 1/256 were prepared in serum-free DMEM. Each dilution was mixed with 100TCID50 of HSV1 and kept for 1h at 37°C in a total volume of 200 µl. Duplicate 100 µl samples were added to HeLa cells in 96-well micro titer plates. DMEM containing 2% heat inactivated newborn-calf serum was added to all inoculated cells including positive and negative controls and were incubated at 37°C. The neutralization titer of each serum

was considered as: the reciprocal titer of the highest dilution preventing CPE in inoculated cells.

# Preparation of recombinant glycoprotein D of HSV1-gD

Plaque purified recombinant baculo virus containing full-length gD-1 gene with titer of 10<sup>7</sup>- 10<sup>8</sup> pfu/ml was applied to infect Spodoptera frugiperda clone 9 (Sf9) insect cells. The cells were grown in Grace's medium containing 10% fetal bovine serum. In order to suitable culture of the get supplementation of the growth medium with 3.5 g/l yeast extract and 3.3 g/l Lactalbumin Hydrolysate was necessary. Gentamycin sulfates (50 mg/ml) and penicillin (100 IU/ml) were also added, and the pH was adjusted to 6.2.

cells were inoculated with purified Sf9 recombinant baculovirus carrying full length of gD-1 gene at multiplicity of infection of 1 and incubated at 27°C for 5 days. The supernatant of the infected cells was collected, filtered through a 0.45 µm pore size filter and stored at -80 °C. In order to produce recombinant gD-1, Sf9 cells were inoculated with purified virus at multiplicity of infection of 5 and incubated at 27 °C for 96h. The cells were harvested and washed three times with ice-cold sterile phosphate buffer (pH 7). The cells were subjected to freezing and thawing for 3 times then were sonicated at 20KHZ, 5 times for 30 seconds(1, 16). The lysate was first treated with CHAPS (0.1mM) and PMSF and then electrophoresed under denaturing conditions and in order to confirm gD-1, Western blotting was done using monoclonal anti HSV1-gD (8).

# Purification of HSV1 by sucrose density gradient centrifugation

Vero cells infected at a multiplicity of infection (moi) of 1 pfu/cell. After appearing the CPEs, the virus was harvested and purified on sucrose density gradients. Briefly, cell debris were removed by centrifugation at 8000g and 4°C for 30 min and the virus was concentrated by centrifugation at 100000 g for 2 hr. The pellet was re-suspended in TNE buffer, pH 7.5 (0.1 M Tris-HCl, 0.1 M NaCl, 0.001 M EDTA) and then purified on step gradient consisting of 30%, 40%, 50% and 60% of sucrose in TNE

and centrifuged at 75000 g for 5 hr. Purified virus was collected, dialyzed and stored at -80°C (16).

### **Immunoblotting**


Virus seeds (either purified virus by sucrose gradient centrifugation recombinant gD-1 protein were mixed separately with electrophoresis sample buffer ( 50mM Tris pH 8.6, 10% glycerol, 2% W/V sodium dodecvl sulfate. and bromophenol blue) and boiled for 10 min. Sodium dodecyl sulfate polyacrylamide gel electrophoresis (PAGE) was performed in 12% polyacrylamide and the separated proteins were transferred to nitrocellulose paper. The proteins bonds on nitrocellulose were reacted with total human HSV1 antibody after blocking in 5% gelatin. Two fold dilutions of positive and negative serum samples for neutralization antibody were

tasted by Western blotting. The bonds were detected by reacting the blots with alkaline phosphate- conjugated anti human antibody and visualized by diaminobenzidine (DAB) as substrate (15).

### **Results**

#### HSV1 specific neutralizing antibodies

Micro neutralization test was performed using

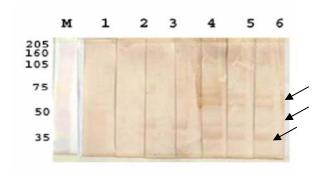

**Fig. 1:** Western blotting results for gD1 purified protein. Its bands are indicated with arrows.

100 TCID50 of the wild-type HSV1. Virus neutralization titer was the reciprocal of the highest dilution of the serum that completely inhibited CPE

in the inoculated cells. Data in Table 1 shows that no HSV1 specific neutralizing antibody (titer < 8) was detected in 11% of the tested sera, while 89% had neutralizing antibody against HSV1.

### **SDS-PAGE** and Western blot analysis

The SDS-PAGE result of the whole virus is showed in Figure 1. The virus specific bands cannot be differentiated from the cellular proteins by SDS-PAGE method. Since some cellular protein have similar molecular weight to viral proteins, therefore Western blotting was used to differentiate these proteins.

Western blot analysis using whole virus as an antigen

**Table 1.** Virus neutralization antibodies to the HSV1 in tested sera.

Virus neutralization antibody titer*	≤2	4	8	16	32	64	128	≥256
Percentage of positive sera **	7%	4%	2%	0%	15%	24%	22%	26%

<sup>\*</sup> Virus neutralization titer was expressed as the reciprocal of the highest dilution that completely inhibited virus CPE in inoculated cells

В

Table 2. The number of positive and negative sera in indicated dilution by VNT and Western blot.

Positive	Negative	Total
7	0	7
1	7	8
8	7	15
	7 1 8	7 0 1 7 8 7

<sup>\*</sup> VNT, virus neutralization test

VNT* Results		
Westown blotting	Positive	N
Western blotting		

Western blotting results	Positive	Negative	Total
Positive	15	1	16
Negative	1	5	6
Total	16	6	22

<sup>\*\*</sup> Two fold dilutions of 100 serum samples were tested against 100 TCID50 of the HSV1.

All of the serum samples with neutralization antibody titers of 32, 64, and 128 were used in Western blotting. Starting with 1000, two fold dilutions were made from the tested sera by VNT and subjected to Western blotting, using whole virus and recombinant HSV-gD1 protein. Since PAGE results showed a non specific band, it was confirmed by Western blotting.

# Western blot analysis using recombinant glycoprotein D of HSV1 (gD-1)

Western blot analysis of baculovirus infected cells showed four distinct bands with apparent molecular weight of ~51,~43.5,~33, and 31 KDa that reacted with monoclonal anti gD-1 antibody (fig.1). The serum samples with positive results in Western blotting using HSV1 were positive at dilutions of 1000, 2000, and 4000 using recombinant gD-1 as antigen (fig. 1).

### Comparison of VNT and Western blotting

In this study, in order to control confounding effects (such as: age, sex ...), the same serum samples were used in both VNT and Western blotting. We performed MC Nemar's test to compare results of VNT and Western blotting (Table 2).

The dilutions of 1000 (A), 4000 (B) of tested sera in WB and 32(A), 128(B) in VN antibody titers were chosen to show their statistical differences, respectively.

Comparing the results as shown in table 2, there is no statistically significant difference.

### **Discussion**

Herpes simplex virus type 1 (HSV1) is a human pathogen which causes a variety of diseases including cold sores, eye and genital infections,

neonatal infections and encephalitis (3). Different methods are used for HSV antibody detection including ELISA, IF, CF and virus neutralization test (VNT). Virus neutralization test is more specific than other serological tests (17). This method requires cell culture and skilled personal and takes long time to perform.

We tried to optimize Western blotting test using HSV1 whole particles and baculovirus derived HSV1-gD.

Based on the results of MC Nemar's statistical test, differences between the two tests are not statistically significant (p>0.05). Western blotting has many advantages compared with VNT and can detect low titers antibody in In this study, 100 serum samples with reciprocal virus neutralization antibody titers of 32, 64, and 128 were chosen and dilutions of 1000, 2000, and 4000 of the sera were prepared and tested by Western blotting, using whole HSV-1 and HSV-gD recombinant proteins. Similarly positive and negative sera were used as controls.

patient's serum. Considering these advantages, Western blotting can be used instead of VNT with similar sensitivity and specifity.

Using recombinant baculo derived HSV1-gD as antigen in Western blotting showed similar results compared with using whole virus. We believe that recombinant proteins should be used instead of whole infectious virus in Western blotting and other serologic tests.

### References

- 1. Handler CG, Eisenberg RJ, Cohen GH. Oligomeric structure of glycoproteins in herpes simplex virus type 1. J Virol. 1996 Sep;70(9):6067-70.
- 2. Wagner EK, Bloom DC. Experimental investigation of herpes simplex virus latency. Clin Microbiol Rev. 1997 Jul;10(3):419-43.
- 3. Cohen GH, Dietzschold B, Ponce de Leon M, Long D, Golub E, Varrichio A, et al. Localization and synthesis of an antigenic determinant of herpes simplex virus glycoprotein D that stimulates the production of neutralizing antibody. J Virol. 1984 Jan;49(1):102-8.
- 4. Baines JD, Roizman B. The UL10 gene of herpes simplex virus type 1 encodes a novel viral glycoprotein, gM, which is present in the virion and in the plasma membrane of infected cells. J of Virol 1993;67:1441-52.
- 5. Hutchinson L, Browne H, Wargent V, Davis-Poynter N, Primorac S, Goldsmith K, et al. A novel herpes simplex virus glycoprotein,

- gL, forms a complex with glycoprotein H (gH) and affects normal folding and surface expression of gH. J Virol. 1992 Apr;66(4):2240-50.
- 6. Burke RL. Development of a herpes simplex virus subunit glycoprotein vaccine for prophylactic and therapeutic use. Rev Infect Dis. 1991 Nov-Dec;13 Suppl 11:S906-11.
- 7. Ghiasi H, Kaiwar R, Nesburn AB, Slanina S, Wechsler SL. Expression of seven herpes simplex virus type 1 glycoproteins (gB, gC, gD, gE, gG, gH, and gI): comparative protection against lethal challenge in mice. J Virol. 1994 Apr;68(4):2118-26.
- 8. Ghiasi H, Nesburn AB, Kaiwar R, Wechsler SL. Immunoselection of recombinant baculoviruses expressing high levels of biologically active herpes simplex virus type 1 glycoprotein D. Arch Virol. 1991;121(1-4):163-78.
- 9. Sisk WP, Bradley JD, Leipold RJ, Stoltzfus AM, Ponce de Leon M, Hilf M, et al. Highlevel expression and purification of secreted forms of herpes simplex virus type 1 glycoprotein gD synthesized by baculovirus-infected insect cells. J Virol. 1994 Feb;68(2):766-75.
- 10. Ishizaka ST, Piacente P, Silva J, Mishkin EM. IgG subtype is correlated with efficiency of passive protection and effector function of anti-herpes simplex virus glycoprotein D monoclonal antibodies. J Infect Dis. 1995 Oct;172(4):1108-11.
- 11. Isola VJ, Eisenberg RJ, Siebert GR, Heilman CJ, Wilcox WC, Cohen GH. Fine

- mapping of antigenic site II of herpes simplex virus glycoprotein D. J Virol. 1989 May;63(5):2325-34.
- 12. Minson AC, Hodgman TC, Digard P, Hancock DC, Bell SE, Buckmaster EA. An analysis of the biological properties of monoclonal antibodies against glycoprotein D of herpes simplex virus and identification of amino acid substitutions that confer resistance to neutralization. J Gen Virol. 1986 Jun;67 (Pt 6):1001-13.
- 13. Muggeridge MI, Wilcox WC, Cohen GH, Eisenberg RJ. Identification of a site on herpes simplex virus type 1 glycoprotein D that is essential for infectivity. J Virol. 1990 Aug;64(8):3617-26.
- 14. Straus SE, Savarese B, Tigges M, Freifeld AG, Krause PR, Margolis DM, et al. Induction and enhancement of immune responses to herpes simplex virus type 2 in humans by use of a recombinant glycoprotein D vaccine. J Infect Dis. 1993 May:167(5):1045-52.
- 15. Soleimanjahi H, Roostaee MH, Rasaee MJ, Mahboudi F, Kazemnejad A, Bamdad T, et al. The effect of DNA priming-protein boosting on enhancing humoral immunity and protecting mice against lethal HSV infections. FEMS Immunol Med Microbiol. 2006;46(1):100-6.
- 16. Deutscher MP. Guide to protein purification. Academic Press San Diego. 1990:665-70.
- 17. Mahy BWJ, Kangro HO. Virology Methods Manual. Academic Press Great Britain, Glasgow 1996:108-9.